


CPS Locator™ 3D Delivery Catheters

DS2C100, DS2C200, DS2C300

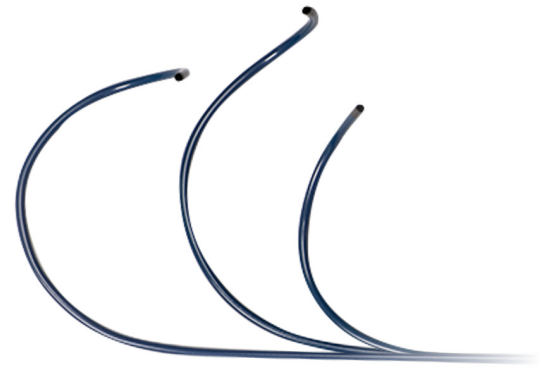

Product Highlights

- The CPS Locator™ 3D Delivery Catheter is a single-use percutaneous catheter intended to introduce pacing or defibrillator leads.
- Each CPS Locator 3D Delivery Catheter offers a unique distal curve to accommodate various anatomies.
- The catheter is equipped with a hemostatic valve, and the distal tip is radiopaque to facilitate imaging under fluoroscopy.
- The Agilis HisPro™ Slitter is designed to ensure alignment and protect the lead during catheter slitting.
- Designed to accommodate a 6F lead with a recommended minimum length of 58 cm for the 42 cm catheters and 65 cm for the 45 cm catheters.¹
- The catheter has an inner diameter of 7 F, an outer diameter of 9 F, and the dilator is compatible with a 0.035" guidewire.

Physical Specifications

MODEL NUMBER	MODEL NAME		OVERALL LENGTH (CM)	WORKING LENGTH (CM)	INNER DIAMETER	OUTER DIAMETER	MINIMUM LEAD LENGTH
DS2C100-42	CPS Locator™ 3D, Small		45 cm	42 cm	7 F (2.44 mm)	9 F (3.00 mm)	58 cm
DS2C200-42	CPS Locator™ 3D, Medium		45 cm	42 cm	7 F (2.44 mm)	9 F (3.00 mm)	58 cm
DS2C300-42	CPS Locator™ 3D, Large		45 cm	42 cm	7 F (2.44 mm)	9 F (3.00 mm)	58 cm
DS2C200-45	CPS Locator™ 3D, Medium, Extra Long		48 cm	45 cm	7 F (2.44 mm)	9 F (3.00 mm)	65 cm
DS2C300-45	CPS Locator™ 3D, Large, Extra Long		48 cm	45 cm	7 F (2.44 mm)	9 F (3.00 mm)	65 cm

Accessories

AVAILABLE SEPARATELY	MODEL NAME	MODEL NUMBER
Slitter	Agilis HisPro™ Slitter	Model DS3A001

Manufactured by:
CenterPoint Systems LLC
3338 Parkway Blvd
West Valley City, UT 84119, U.S.A.
Tel: +1-877-848-0828
Info@centerpoint-systems.com

Distributed by:
Abbott Medical
15900 Valley View Court
Sylmar, CA 91342 USA
+1-818-362-6822

U.S. Customer Support:
1-800-722-3774

References

1. Abbott. Data on File, Technical Report. Item: 91112329

Brief Summary: Prior to using these devices, please review the Instructions for Use for a complete listing of indications, contraindications, warnings, precautions, potential adverse events and directions for use.

[™] Indicates a trademark of the Abbott group of companies.

[®] Indicates a third-party trademark, which is property of its respective owner.

Abbott
15900 Valley View Court, Sylmar, CA 91342
Tel: +1 818 362 6822
Abbott.com

© 2025 Abbott. All Rights Reserved.
MAT-2503989 v1.0 | Item approved for OUS use.

