


# CPS Locator™ 3D Catheters






DS2C100, DS2C200, DS2C300

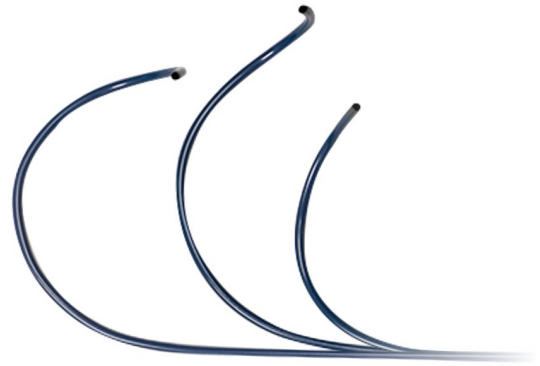

## Product Highlights

- The CPS Locator™ 3D is a single-use percutaneous catheter intended to introduce various types of catheters and pacing or defibrillator leads.
- Each CPS Locator™ 3D catheter offers a unique distal curve to accommodate various anatomies.
- The catheter is equipped with a hemostatic valve, and the distal tip is radiopaque to facilitate imaging under fluoroscopy.
- Slitter for CPS Locator™ 3D catheter is designed to ensure alignment and protect the lead during catheter slitting.
- Accommodates a 6 F lead with a minimum length of 58 cm for the 42 cm catheters and 65 cm for the 45 cm catheters.
- The catheter has an inner diameter of 7 F, an outer diameter of 9 F, and the dilator is compatible with a 0.035" guidewire.

## Physical Specifications

MODEL NUMBER	MODEL NAME		OVERALL LENGTH (CM)	WORKING LENGTH (CM)	INNER DIAMETER	OUTER DIAMETER	MINIMUM LEAD LENGTH
DS2C100-42	CPS Locator™ 3D, Small		45 cm	42 cm	7 F (2.44 mm)	9 F (3.00 mm)	58 cm
DS2C200-42	CPS Locator™ 3D, Medium		45 cm	42 cm	7 F (2.44 mm)	9 F (3.00 mm)	58 cm
DS2C300-42	CPS Locator™ 3D, Large		45 cm	42 cm	7 F (2.44 mm)	9 F (3.00 mm)	58 cm
DS2C200-45*	CPS Locator™ 3D, Medium, Extra Long		48 cm	45 cm	7 F (2.44 mm)	9 F (3.00 mm)	65 cm
DS2C300-45*	CPS Locator™ 3D, Large, Extra Long		48 cm	45 cm	7 F (2.44 mm)	9 F (3.00 mm)	65 cm

\*Please check with your local Abbott representative for product availability

## Accessories

INCLUDED	AVAILABLE SEPARATELY
Slitter	Agilis HisPro™ Slitter (Model DS3A001)



**Manufactured by:**

CenterPoint Systems LLC  
3338 Parkway Blvd  
West Valley City, UT 84119, U.S.A.  
Tel: +1-877-848-0828  
Info@centerpoint-systems.com

**Distributed by:**

Abbott Medical  
15900 Valley View Court  
Sylmar, CA 91342 USA  
+1-818-362-6822

**U.S. Customer Support:**

1-800-722-3774

**Rx Only**

**Brief Summary:** Prior to using these devices, please review the Instructions for Use for a complete listing of indications, contraindications, warnings, precautions, potential adverse events and directions for use.

™ Indicates a trademark of the Abbott group of companies.

‡ Indicates a third party trademark, which is property of its respective owner.

© 2023 Abbott. All Rights Reserved.

MAT-2301332 v3.0 | Item approved for U.S. use only.

