



AMG PP Adult Oxygenator

PRODUCT HIGHLIGHTS

Microporous hollow-fiber oxygenator with an integral heat exchanger for use in extracorporeal circuits.

Key Features

- Polypropylene (PP) fibers
- Phosphorylcholine (PC) coating
- Integrated arterial and venous luer ports
- Compatible with any circuit with 3/8" tubing
- Includes two pigtailed stopcocks and a gas line with filter
- Compatible with the same holder used with the AMG PMP Adult Oxygenator

PRODUCT SPECIFICATIONS

PART NUMBER	DESCRIPTION
US5400	AMG PP Adult Oxygenator Box of 4



AMG PP
Adult Oxygenator


TECHNICAL SPECIFICATIONS

TECHNICAL SPECIFICATION	AMG PP ADULT OXYGENATOR
Maximum Blood Flow Rate	7.0 L/min
Priming Volume	220 mL
Gas Exchange Membrane Surface Area	1.81 m ²
Heat Exchanger Surface Area	0.08 m ²
Gas Exchange Material (Membrane Type)	Polypropylene (PP)
Heat Exchange Material	Stainless Steel
Surface Coating	Phosphorylcholine (PC)
Inlet and Outlet	3/8" Barbed Connectors

PRESSURE, HEAT EXCHANGE AND GAS EXCHANGE CURVES


O₂ TRANSFER

Efficiency of blood oxygenation

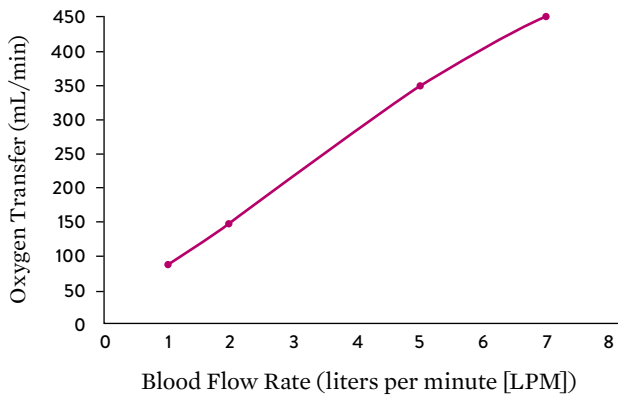


CO₂ TRANSFER

Efficiency of blood carbon dioxide removal

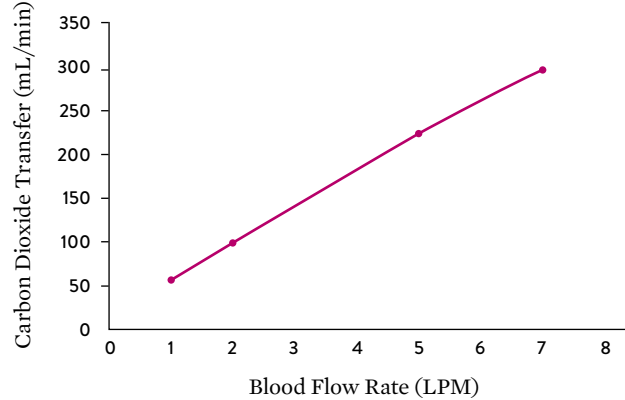


PRESSURE DROP

Difference in blood pressure before and after the oxygenator

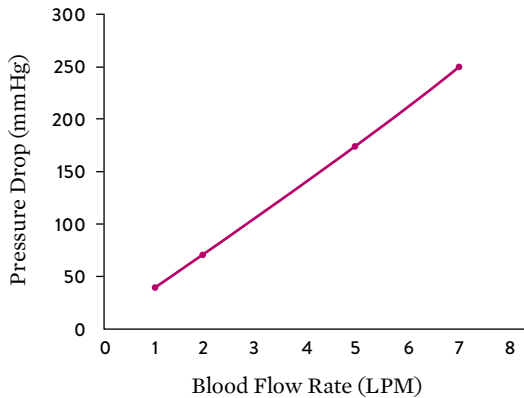

HEAT EXCHANGE

Efficiency of heat exchange

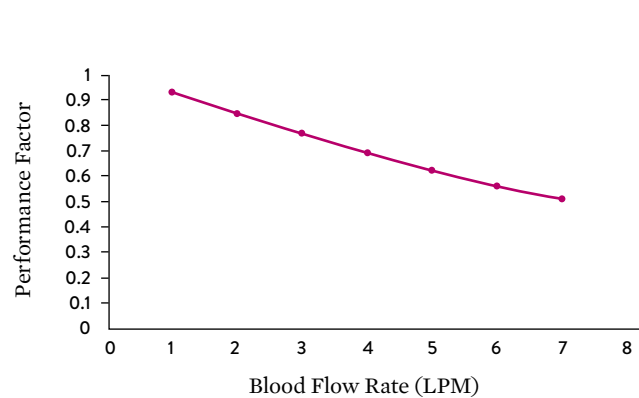

AMG PP Adult Oxygenator and Oxygenator Holder are manufactured by Eurosets s.r.l. and distributed by Abbott.

Abbott

6101 Stoneridge Dr., Pleasanton, CA 94588 USA, Tel: 1 925 847 8600
Cardiovascular.Abbott/CentriMag

Rx Only

Brief Summary: Prior to using these devices, please review the Instructions for Use for a complete listing of indications, contraindications, warnings, precautions, potential adverse events and directions for use.

AMG PP Adult Oxygenator Temporary Expanded Indication: A.M.G. MODULE PP NO T.P. STERILE is indicated for use as part of an ECMO circuit for longer than 6 hours to treat patients with acute respiratory failure and/or acute cardiopulmonary failure during the COVID-19 public health emergency.

AMG PP Adult Oxygenator Indications [510(k) Clearance; 6-hour use]: The A.M.G. MODULE PP NO T.P. STERILE is intended for use in adult surgical procedures requiring extracorporeal gas exchange support and blood temperature control for periods of up to 6 hours. The A.M.G. MODULE PP NO T.P. STERILE is intended to be used in an extracorporeal perfusion circuit to oxygenate and remove carbon dioxide from the blood and to cool or warm

the blood during routine cardiopulmonary bypass procedures up to 6 hours in duration. Contact with blood for a longer period of time is inadvisable. A.M.G. MODULE PP NO T.P. STERILE is an adult oxygenator. The A.M.G. MODULE PP NO T.P. STERILE should be used in combination with medical devices listed in section "Medical devices for use with the A.M.G. MODULE PP NO T.P. STERILE."

AMG PP Adult Oxygenator Contraindications [510(k) Clearance; 6-hour use]: Do not use the device for any purpose other than indicated. This device used for any other purposes than for the indicated intended use is the responsibility of the user. Relative contraindications to CPB procedures applies.

[™] Indicates a trademark of the Abbott group of companies.
‡ Indicates a third party trademark, which is property of its respective owner.

© 2022 Abbott. All Rights Reserved.
MAT-2004391 v2.0 | Item approved for U.S. use only

